
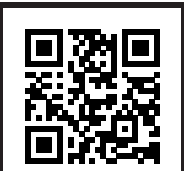
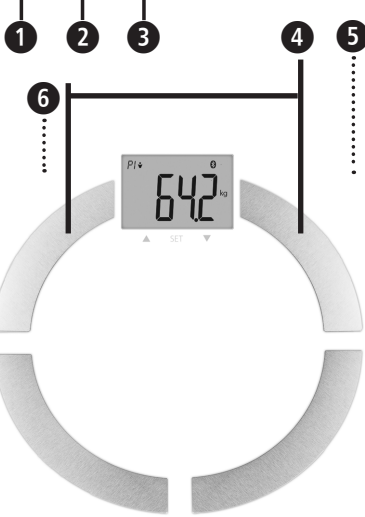
medisana®

BS 444

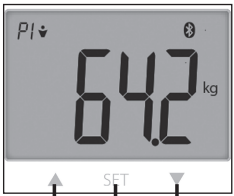Manual

Service

medisana GmbH
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

REF 40444

EAN 40 15588 40444 3

Körperanalysewaage BS 444

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Body analysis scale BS 444

- ▲ Up button
- SET button
- ▼ Down button
- Electrodes
- Battery compartment (on the underside)
- Switch for weight units (on the underside) for kilograms (kg), pounds (lb), stones (st)

Man

Man (athletic)

Woman

Woman (athletic)

Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://docs.medisana.com/40444>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!
 - This scale is not suitable for people with pacemakers or other medical implants. The indication of body fat may be inaccurate in the event of diabetes or other medical/physical conditions. The same applies to highly trained sportsmen.
 - This device is not suitable for pregnant women!
 - Any treatment or diet for overweight or underweight requires qualified advice from appropriate specialists (doctor, dietician). The values determined with the scale can be of help.
 - Risk of tilting! Never place the scales on an uneven surface. Never step onto just one corner of the scale. While weighing, never stand on one side or on the edge of the scale.
 - Danger of slipping! Never step onto the scale with wet feet. Never step onto the scale while wearing socks.
 - Use the device only for its intended purpose, as directed by the user instructions. In the event of misuse, all guarantees become null and void.
 - The scale is manufactured for home use. It is not intended for commercial use in hospitals or other medical facilities.
 - Do not use the device if it does not work properly, if it has been dropped, dropped into water or damaged.
 - The scale has a measuring range up to 180 kg / 396 lb / 28 st 4 lb. Do not overload the scale.
 - Place the scale on a firm level surface. Soft, uneven surfaces are unsuitable for measurement and lead to incorrect results.
 - Place the scale in a location where neither extreme temperatures nor extreme humidity are present.
 - Keep the scale away from water.
 - Please handle the scale with care. Avoid shocks and vibrations of the scale. Do not drop it.
 - Stand carefully on the scale. Do not hop or jump on the tread surface, as this may damage the weighing mechanism.
 - Do not attempt to disassemble the device, otherwise the warranty will become void. This device does not contain any parts that can be serviced or replaced by the user.
 - In the event of malfunctions, do not proceed with any unauthorised repair of the device. Otherwise, this will render any warranty claims void. Have repairs performed only by an authorised service centre.
- Battery safety instructions**
- Do not take the batteries apart!
 - If required, clean the battery and device contacts before insertion!
 - Removed spent batteries from the device immediately!
 - Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! In the event of contact with battery acid, immediately rinse the affected area with plenty of clean water and consult a doctor immediately!
 - If a battery is swallowed, consult a doctor immediately!

- Insert the batteries correctly, paying attention to the polarity!
- Keep the battery compartment closed properly!
- If the device is not used for an extended period of time, remove the batteries!
- Keep the batteries out of the reach of children!
- Do not re-charge the batteries! There is a risk of explosion!
- Do not short-circuit the device! There is a risk of explosion!
- Do not throw into a fire! There is a risk of explosion!
- Store unused batteries in the packaging and not near metal objects to avoid a short circuit!

Simple weighing

- Stand on the scale and stand still. The device automatically switches itself on.
- Your weight is measured, flashes twice and then appears constantly.
- Step down off the scale.

Optional: Transfer the measured weight* to the VitaDock+ App.

Prerequisites:

- Bluetooth is activated on your Android or iOS device
 - The VitaDock+ App has been started.
 - The scale is already listed in the device list in the App.
- * No body analysis data are transmitted here and no weight history is created.

Body analysis

The scale can store the data of 8 persons: gender, age, height.

- Briefly press the middle of the scale surface with your foot. The display will show "0.0 kg".
- Long press the SET button ➊. The memory location flashes in the display.
- Now you can select the desired memory location for your user profile (0 - 8) by pressing button ➊ or button ➋.
- NOTE: If you use the scale in guest mode "P0", there will be no data storage. Transmission of the measured values via Bluetooth is not possible in this mode.
- Press the SET button ➋ to save the setting.
- Continue to adjust the gender, athlete mode, size, and age settings.
- After you have set your personal data, the display shows "0.0". To finish the scale settings, simply stand barefoot on the scale.
- The display shows the following measured values:
 - Weight • Proportion of body fat • Proportion of body water • Proportion of muscle mass • Bone weight • BMR (calorie requirement) • BMI

When a complete measurement could be taken, the data are then sent to the VitaDock+ App via Bluetooth.

Guidelines and standards

We hereby declare that this device is in compliance with the essential requirements of the European RED 2014/53/EU directive. The complete EU declaration of conformity can be requested from medisana GmbH, Carl-Schurz-Str. 2, 41460 Neuss, Germany.

Voltage supply: 3 batteries (Type AAA, 1.5 V, Micro LR03)

Disposal: This product may not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

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